



We're here to help.

Your Checklist for Preparing to Order the Correct SleepSafe® Bed

When ordering your *SleepSafe® Bed*, please select the following choices:

- 1. Rail Height:** LOW MEDIUM TALL TALL EXT Open Closed
- 2. Foundation:** FIXED Box Spring Bunkie Board
 ELECTRIC HiLo Combo Articulation
 MANUAL HiLo HiLo Combo Articulation Electric Articulation / Manual HiLo
- 3. Finish:** OAK MAPLE CHERRY MAHOGANY ALDER
 SOLID COLOR White Blue Green Red Yellow Orange Purple Pink Lt. Brown Dk. Brown Gray Black
 MULTI COLOR (Please call for a form to select your color combination.)
- 4. Bed Size:** TWIN FULL QUEEN
- 5. Padding:** PADDING AROUND WINDOWS PADDING OVER WINDOWS NO PADDING
PADDING COLOR Cream Gray Light Blue Taupe Dark Blue Pink
- 6. Windows:** CLEAR WINDOWS MESH WINDOWS CLEAR/MESH COMBINATION
- 7. Extras:** HEADBOARD WINDOW FOOTBOARD WINDOW JENX DREAMA SYSTEM
 IV POLE KIT TUBING CHANNEL ACCESS DUAL VIEW (Tall Bed)

All *SleepSafe® Beds* include these items:

- Bed Enclosure
- Premium Mattress
- Assembly Instructions
- Standard Shipping Within the Continental US
- Foundation
- Casters
- Tools for Assembly



3629 Reed Creek Drive · Bassett, VA 24055 · Local: 276-627-0088 · Fax: 276-627-0234



SleepSafeBed.com | 866.852.2337

Safety First

SleepSafe® Beds are a proven, ultra-safe solution to address problems with entrapment and falls. Unique in the industry, these beds meet or exceed FDA guidelines on the **7 Zones of Entrapment**.

THREE MODELS of SleepSafe® Beds are available in TWIN, FULL and QUEEN sizes with safety rail height protection ranging from 8¾" to 53" above the mattress depending on the configuration. Appropriate SAFETY RAIL HEIGHT is based on the activity level of the user, with more active users requiring additional height.




*SleepSafe® II
Medium Bed
with Alder finish*

Easily Moved

Casters provide mobility and can be locked.

Premium Memory Foam Mattress

SleepSafe® Beds use CertiPUR-US® compliant foam products for our mattresses featuring gel technology. They provide optimized pressure redistribution, rapid heat dissipation, and increases in breathability compared to standard foam. Mattress coverings are fluid resistant and support up to 350 pounds.

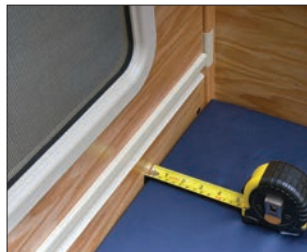

Resists Burrowing

The **SleepSafe®** precision-fit custom mattress fits tightly within the bed frame. Even with aggressive compression, it is difficult to burrow deeply enough to create an entrapment issue or safety hazard.

Safest Yet!

SleepSafer® Extension Conversion Kit increases bed safety height another 16¾". The Kit, with or without top enclosure, may be included with a new bed order or can be sent as a retrofit for most existing beds. Screen mesh windows offer additional ventilation.

Easy-access Safety Rails

Full-length wood safety rails are both hinged and removable for quick and easy accessibility.

Choose Clear or Mesh Windows

Clear, break-resistant windows, or our new mesh windows, provide visibility for both the occupant and caregiver. The need to see and be seen is important for comfort and safety.

TRANSFER HEIGHTS

Please consider the **TRANSFER HEIGHT** before ordering your bed. It is the distance from the top of the side panel (below the safety rail) to the floor. SleepSafe® - LOW - 20¾" or 22¼" Transfer Height
SleepSafe® II - MEDIUM - 30½" or 32" Transfer Height
SleepSafer® - TALL - 24½" or 26" Transfer Height




SleepSafeBed.com | 866.852.2337

<i>SleepSafe® Bed Specifications</i>	<i>Low Bed SleepSafe®</i>	<i>Medium Bed SleepSafe® II</i>	<i>Tall Bed / Extension Bed SleepSafer®</i>
FOUNDATION SAFETY RAIL HEIGHTS <i>The distance from the top of mattress to the top of the safety rail.</i>			
Fixed w/Box Spring (Frame attachment choice)	12" or 8¾"	17¾", 14" or 10"	31" or 24¾" / 47¾" or 41½"
Fixed w/Bunkie Board (Frame attachment choice)	17¼" or 14"	27" or 21¾"	36¼" / 53"
Manual Articulation (Frame attachment choice)	15" or 11¾"	21", 17" or 13¼"	34" or 27¾" / 50¾" or 44½"
Manual Electric Combo (LO to HI Settings based on frame attachment)	LO 11½" to EVEN with rail	LO 21¼" to HI 9½" or LO 16" to HI 4¼"	LO 30½" to HI 18¾" / LO 47½" to HI 35¾"
Manual Hi-Lo (LO to HI Settings based on frame attachment)	LO 11½" to EVEN with rail	LO 21¼" to HI 9½" or LO 16" to HI 4¼"	LO 30½" to HI 18¾" / LO 47½" to HI 35¾"
Manual Hi-Lo Combo (LO to HI Settings based on frame attachment)	LO 11½" to EVEN with rail	LO 21¼" to HI 9½" or LO 16" to HI 4¼"	LO 30½" to HI 18¾" / LO 47½" to HI 35¾"
Electric Articulation (Frame attachment choice)	15" or 11¾"	21", 17" or 13¼"	34" or 27¾" / 50¾" or 44½"
Electric Hi-Lo Combo (LO to HI Settings based on frame attachment)	LO 11½" to HI 1¼" below rail	LO 21¼" to HI 11" or LO 16" to HI 5¾"	LO 30½" to HI 20¾" / LO 47½" to HI 37¼"
FOUNDATION MATTRESS TO FLOOR HEIGHTS <i>The distance from the top of the mattress to the floor.</i>			
Fixed w/Box Spring (Frame attachment choice)	23¾" or 27"	28", 32" or 35½"	23¾" or 30"
Fixed w/Bunkie Board (Frame attachment choice)	18½" or 21¾"	18¾" or 24"	18½"
Manual Articulation (Frame attachment choice)	20¾" or 24"	25", 28¾" or 32½"	20¾" or 27"
Manual Electric Combo (LO to HI Settings based on frame attachment)	LO 24¼" to HI 36"	LO 24½" to HI 36¼" or LO 29¾" to HI 41½"	LO 24¼" to HI 36"
Manual Hi-Lo (LO to HI Settings based on frame attachment)	LO 24¼" to HI 36"	LO 24½" to HI 36¼" or LO 29¾" to HI 41½"	LO 24¼" to HI 36"
Manual Hi-Lo Combo (LO to HI Settings based on frame attachment)	LO 24¼" to HI 36"	LO 24½" to HI 36¼" or LO 29¾" to HI 41½"	LO 24¼" to HI 36"
Electric Articulation (Frame attachment choice)	20¾" or 24"	25", 28¾" or 32½"	20¾" or 27"
Electric Hi-Lo Combo (LO to HI Settings based on frame attachment)	LO 24¼" to HI 34½"	LO 24½" to HI 34¾" or LO 29¾" to HI 40"	LO 24¼" to HI 34½"
FLOOR TO RAIL HEIGHT <i>The distance from the floor to the top of the rail.</i>	35¾"	45¾"	54¾" / 72¾"
OVERALL OUTSIDE DIMENSIONS			
Twin	42¼" W x 80½" L x 40¾" H	42¼" W x 80½" L x 50¾" H	42¼" W x 80¼" L x 56" H (72¾" H w/Ext.)
Full	56" W x 80" L x 40¾" H	56" W x 80" L x 50¾" H	56" W x 80" L x 57" H (72¾" H w/Ext.)
Queen	65" W x 80" L x 40¾" H	65" W x 80" L x 50¾" H	65" W x 80" L x 57" H (72¾" H w/Ext.)
MATTRESS SIZE <i>Width without pads / Width with pads / Length</i>			
Twin / Full / Queen	38" W / 36" W / 77" L	52" W / 50" W / 77" L	61" W / 60" W / 77" L
UNDER BED FLOOR CLEARANCE			
Standard Models	8"	8"	7¾"
TRANSFER HEIGHT <i>The distance from top edge of the folded down safety rail to the floor.</i>	22¼" with Safety Rail On 20¾" with Safety Rail Off	32" with Safety Rail On 30½" with Safety Rail Off	26" with Safety Rail On 24½" with Safety Rail Off



Photo Reference Guide for Bed Measurements

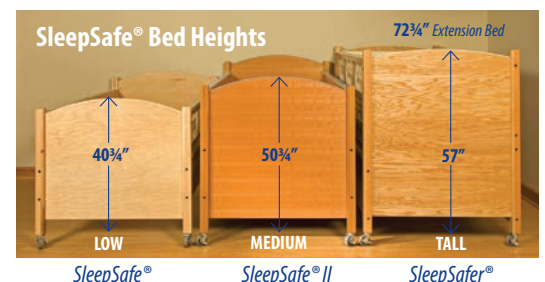



*SleepSafe II
Medium Bed
with Mahogany Finish*

Foundations

FIXED – The mattress rests on a fixed foundation, which can be a box spring or bunkie board, depending on the safety rail protection needed. A bunkie board allows 6" of additional safety rail height.



ARTICULATED – The mattress rests on an ELECTRIC or MANUAL frame that has both head and knee articulation, operated by a hand-held remote or crank.

HI-LO – The mattress rests on an ELECTRIC frame that offers mattress height adjustability as well as head and knee articulation. The surface of the mattress can be adjusted from 22¾" to 38¾" from the floor.



Break-resistant windows are also available for the headboard and/or footboard.

Finish Choices

SleepSafe® Beds can be finished in wood or color choices:



SOLID COLORS / MULTICOLOR

White, Blue, Green, Red, Yellow, Orange, Purple, Pink, Black, Dk. Brown, Lt. Brown and Gray



*SleepSafe II
Medium Bed
with Oak Finish*

Accessories

DUAL VIEW SAFETY RAILS

For more visibility or accessibility, our SleepSafer® TALL bed may be ordered with see-through rails on **both** sides.

DREAMA® POSITIONING SYSTEM

SleepSafe® Beds can provide a retrofit for the DREAMA™ Positioning System by Jenx, a support system for effective 24-hour postural management, suitable for children and adults.

PADDING

Custom-fit padding provides protection from contact with the bed frame. The ¾" high-density foam is covered with easy to clean waterproof nylon. Pads can be ordered to surround the windows or cover the windows. They can be removed for cleaning.

IMPORTANT: When required, padding needs to be indicated at the time of ordering your bed.

IV ATTACHMENT KIT

A stainless steel IV pole, along with an access area, may be ordered with your bed.



SleepSafeBed.com | 866.852.2337