

We're here to help.



STELLA, in a picture shared by her mom Sarah, is one of our SleepSafe SuperStars. Discover more at #SleepSafeSuperStar!

Your checklist for preparing to order the correct SleepSafe® Bed

When ordering your **SleepSafe® Bed**, please select the following choices:

- 1. Rail Height:** ☐ LOW ☐ MEDIUM ☐ TALL ☐ EXTENSION BED ☐ Open ☐ Closed
- 2. Foundation:** ☐ FIXED ☐ Box Spring ☐ Bunkie Board
☐ ELECTRIC ☐ HiLo Combo ☐ Articulation
☐ MANUAL ☐ HiLo ☐ HiLo Combo ☐ Articulation ☐ Electric Articulation / Manual HiLo
- 3. Finish:** ☐ OAK ☐ MAPLE ☐ CHERRY ☐ MAHOGANY ☐ ALDER
☐ SOLID COLOR ☐ White ☐ Blue ☐ Green ☐ Red ☐ Yellow ☐ Orange ☐ Purple ☐ Pink ☐ Lt. Brown ☐ Dk. Brown ☐ Gray ☐ Black
☐ MULTI COLOR (Please see our website or call for a form to select your color combination.)
- 4. Bed Size:** ☐ TWIN ☐ FULL ☐ QUEEN
- 5. Padding:** ☐ PADDING AROUND WINDOWS ☐ PADDING OVER WINDOWS ☐ NO PADDING
PADDING COLOR ☐ Cream ☐ Gray ☐ Light Blue ☐ Taupe ☐ Dark Blue ☐ Pink
- 6. Windows:** ☐ CLEAR WINDOWS ☐ MESH WINDOWS ☐ CLEAR/MESH COMBINATION
- 7. Extras:** ☐ HEADBOARD WINDOW ☐ FOOTBOARD WINDOW ☐ DUAL VIEW (Tall Bed)
☐ IV POLE KIT ☐ TUBING CHANNEL ACCESS ☐ CUSTOM BED LENGTH
☐ JENX DREAMA POSITIONING SYSTEM ☐ LOW AIR LOSS MATTRESS

All SleepSafe® Beds include:

- Bed Enclosure
- Premium Mattress
- Assembly Instructions and Tools
- Standard Shipping Within Continental US
- Foundation
- Casters



3629 Reed Creek Drive · Bassett, VA 24055 · Local: 276-627-0088 · Fax: 276-627-0234



SleepSafeBed.com | 866.852.2337

Safety First

SleepSafe® Beds are a proven, ultra-safe solution to address problems with entrapment and falls. Unique in the industry, these beds meet or exceed FDA guidelines on the **7 Zones of Entrapment**.

FOUR MODELS of SleepSafe® Beds are available in TWIN, FULL and QUEEN sizes with safety rail height protection ranging from 8¾" to 53" above the mattress depending on the configuration. Appropriate SAFETY RAIL HEIGHT is based on the activity level of the user, with more active users requiring additional height.



*SleepSafe® II
Medium Bed
with Alder finish*

Easily Moved


Casters provide mobility and can be locked.

Premium Memory Foam Mattress

SleepSafe® Beds use CertiPUR-US® compliant foam products for our mattresses featuring gel technology. They provide optimized pressure redistribution, rapid heat dissipation, and increases in breathability compared to standard foam. Mattress coverings are fluid resistant and support up to 350 pounds.

New!

SuperScreen Top enclosure is now standard for Extension Beds.

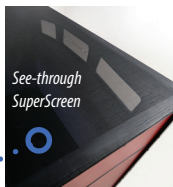

*See-through
SuperScreen*



Safest Yet!


The **SleepSafe® Extension Bed** increases bed safety height another 16¾" above the safety rail height of the SleepSafer® Bed. The bed may be ordered with or without a top enclosure. A retrofit kit for most existing SleepSafer® Beds is available.

Easy-access Safety Rails

Full-length wood safety rails are both hinged and removable for quick and easy accessibility.

Choose Clear or Mesh Windows

Clear, break-resistant windows, or our new mesh windows, provide visibility for both the occupant and caregiver. The need to see and be seen is important for comfort and safety.

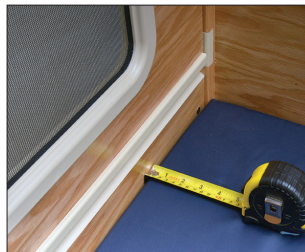

Resists Burrowing

The **SleepSafe®** precision-fit custom mattress fits tightly within the bed frame. Even with aggressive compression, it is difficult to burrow deeply enough to create an entrapment issue or safety hazard.

TRANSFER HEIGHTS

Please consider the **TRANSFER HEIGHT** before ordering your bed. It is the distance from the top of the side panel (below the safety rail) to the floor. SleepSafe® - LOW - 20¾" or 22¼" Transfer Height
SleepSafe® II - MEDIUM - 30½" or 32" Transfer Height
SleepSafer® - TALL - 24½" or 26" Transfer Height



SleepSafeBed.com | 866.852.2337

SleepSafe® Bed Model Measurements

Use this chart to compare the different measurement specifications from model to model.


SleepSafe Bed, LLC
3629 Reed Creek Dr.
Bassett, VA 24055

Toll Free: 866.852.2337
SleepSafeBed.com

<i>SleepSafe® Bed</i> Specifications	<i>Low Bed</i> SleepSafe®	<i>Medium Bed</i> SleepSafe® II	<i>Tall Bed / Extension Bed</i> SleepSafer®
FOUNDATION SAFETY RAIL HEIGHTS <i>The distance from the top of mattress to the top of the safety rail.</i>			
Fixed w/Box Spring (Frame attachment choice)	12" or 8¾"	17¾", 14" or 10"	31" or 24¾" / 47¾" or 41½"
Fixed w/Bunkie Board (Frame attachment choice)	17¼" or 14"	27" or 21¾"	36¼" or 30" / 53" or 46¾"
Manual Articulation (Frame attachment choice)	15" or 11¾"	21", 17" or 13¼"	34" or 27¾" / 50¾" or 44½"
Manual Electric Combo (LO to HI Settings based on frame attachment)	LO 11½" to EVEN with rail	LO 21¼" to HI 9½" or LO 16" to HI 4¼"	LO 30½" to HI 18¾" / LO 47½" to HI 35¾"
Manual Hi-Lo (LO to HI Settings based on frame attachment)	LO 11½" to EVEN with rail	LO 21¼" to HI 9½" or LO 16" to HI 4¼"	LO 30½" to HI 18¾" / LO 47½" to HI 35¾"
Manual Hi-Lo Combo (LO to HI Settings based on frame attachment)	LO 11½" to EVEN with rail	LO 21¼" to HI 9½" or LO 16" to HI 4¼"	LO 30½" to HI 18¾" / LO 47½" to HI 35¾"
Electric Articulation (Frame attachment choice)	15" or 11¾"	21", 17" or 13¼"	34" or 27¾" / 50¾" or 44½"
Electric Hi-Lo Combo (LO to HI Settings based on frame attachment)	LO 11½" to HI 1¼" below rail	LO 21¼" to HI 11" or LO 16" to HI 5¾"	LO 30½" to HI 20¾" / LO 47½" to HI 37¼"
FOUNDATION MATTRESS TO FLOOR HEIGHTS <i>The distance from the top of the mattress to the floor.</i>			
Fixed w/Box Spring (Frame attachment choice)	23¾" or 27"	28", 32" or 35½"	23¾" or 30"
Fixed w/Bunkie Board (Frame attachment choice)	18½" or 21¾"	18¾" or 24"	18½" or 24¾"
Manual Articulation (Frame attachment choice)	20¾" or 24"	25", 28¾" or 32½"	20¾" or 27"
Manual Electric Combo (LO to HI Settings based on frame attachment)	LO 24¼" to HI 36"	LO 24½" to HI 36¼" or LO 29¾" to HI 41½"	LO 24¼" to HI 36"
Manual Hi-Lo (LO to HI Settings based on frame attachment)	LO 24¼" to HI 36"	LO 24½" to HI 36¼" or LO 29¾" to HI 41½"	LO 24¼" to HI 36"
Manual Hi-Lo Combo (LO to HI Settings based on frame attachment)	LO 24¼" to HI 36"	LO 24½" to HI 36¼" or LO 29¾" to HI 41½"	LO 24¼" to HI 36"
Electric Articulation (Frame attachment choice)	20¾" or 24"	25", 28¾" or 32½"	20¾" or 27"
Electric Hi-Lo Combo (LO to HI Settings based on frame attachment)	LO 24¼" to HI 34½"	LO 24½" to HI 34¾" or LO 29¾" to HI 40"	LO 24¼" to HI 34½"
FLOOR TO RAIL HEIGHT <i>The distance from the floor to the top of the rail.</i>	35¾"	45¾"	54¾" / 72"
OVERALL OUTSIDE DIMENSIONS			
Twin	42¼" W x 80½" L x 40¾" H	42¼" W x 80½" L x 50¾" H	42¼" W x 80¼" L x 57" H (72¾" H w/Ext.)
Full	56" W x 80" L x 40¾" H	56" W x 80" L x 50¾" H	56" W x 80" L x 57" H (72¾" H w/Ext.)
Queen	65" W x 80" L x 40¾" H	65" W x 80" L x 50¾" H	65" W x 80" L x 57" H (72¾" H w/Ext.)
MATTRESS SIZE <i>Width without pads / Width with pads / Length</i>			
Twin / Full / Queen	38" W / 36" W / 77" L	52" W / 50" W / 77" L	61" W / 60" W / 77" L
UNDER BED FLOOR CLEARANCE			
Standard Models	8"	8"	7¾"
TRANSFER HEIGHT <i>The distance from top edge of the folded down safety rail to the floor.</i>	22¼" with Safety Rail On 20¾" with Safety Rail Off	32" with Safety Rail On 30½" with Safety Rail Off	26" with Safety Rail On 24½" with Safety Rail Off



Photo Reference Guide for Bed Measurements

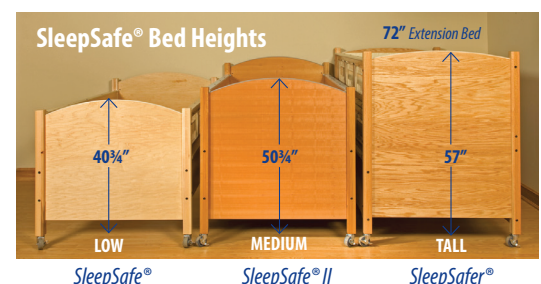

Accessories

PADDING

Custom-fit padding provides a layer of protection from the bed frame. Made with 3/4" high-density foam, it is covered with water-resistant nylon. Pads can either surround the windows or cover the windows.

IMPORTANT: When required, padding needs to be indicated when placing the order for your bed.

DREAMA® POSITIONING SYSTEM

A 24-hour postural management support system that can be used in a SleepSafe bed, and in most cases, be retrofitted to an existing bed.

LOW AIR LOSS MATTRESS

An Alternating Air-pressure Mattress from Blue Chip Medical, can now be fitted for use in a SleepSafe® Bed.

DUAL VIEW SAFETY RAILS

May be ordered for a SleepSafer Tall Bed with safety rails on **both** sides for dual access and visibility.

IV KIT / TUBING CHANNEL

An IV pole kit and a 2"x18" channel for medical tubing, may be ordered with your bed.



*SleepSafe II
Medium Bed
with Mahogany Finish*

Foundations

FIXED – The mattress rests on a fixed foundation, which can be a box spring or bunkie board, depending on the safety rail protection needed. A bunkie board allows 6" of additional safety rail height.

ARTICULATED – The mattress rests on an ELECTRIC or MANUAL frame that has both head and knee articulation, operated by a hand-held remote or crank.

HI-LO – The mattress rests on an ELECTRIC frame that offers mattress height adjustability as well as head and knee articulation. The surface of the mattress can be adjusted from 22¾" to 38¾" from the floor.



Break-resistant windows are also available for the headboard and/or footboard.

Finish Choices

SleepSafe® Beds can be finished in wood or color choices:



SOLID COLORS / MULTICOLOR

White, Blue, Green, Red, Yellow, Orange, Purple, Pink, Black, Dk. Brown, Lt. Brown and Gray



*SleepSafe II
Medium Bed
with Oak Finish
Cream Padding*

PADDING COLORS

