

Jon-Blaise with his big sister Eloise, is one of our SleepSafe SuperStars. Discover more at #SleepSafeSuperStar!

	dist for preparing to order the correct SleepSaf	fe [®] Bed		
When ordering y	our SleepSafe® Bed, please select the following choices:	All SleepSafe® Beds include:		
1. Rail Height:	○ LOW ○ MEDIUM ○ TALL ○ EXTENSION BED ○ Open ○ Closed ○ Screen ○ Wood	Bed Enclosure Premium Mattress Casters		
2. Foundation:	○ FIXEDBox SpringBunkie Board	Assembly Instructions and Tools Standard Shipping Within Continental US		
	© ELECTRICHiLo ComboArticulation Only	11. 3		
	MANUAL HiLo Only HiLo Combo Articulation Only Electric Articulation	& Manual HiLo		
3. Finish:	○ OAK ○ MAPLE ○ CHERRY ○ MAHOGANY ○ ALDER			
	SOLID COLOR _White _Blue _Green _Red _Yellow _Orange _Purple _Pink _Lt. Brov	vn _Dk. Brown _Gray _Black		
	○ MULTI COLOR (<i>Please see our website or call for a form to select your color combination.</i>)			
4. Bed Size:	○TWIN ○ FULL ○ QUEEN			
5. Padding:	 ○ PADDING AROUND WINDOWS ○ PADDING OVER WINDOWS ○ NO PADDING PADDING COLOR Cream Gray Light Blue Taupe Dark Blue Pink 			
6. Windows:	○ CLEAR WINDOWS ○ MESH WINDOWS ○ CLEAR/MESH COMBINATION ○ NO WINDOW	WS		
7. Extras:	 ○ HEADBOARD WINDOW ○ FOOTBOARD WINDOW ○ DUAL VIEW (Tall Bed) ○ IV POLE KIT ○ TUBING CHANNEL ACCESS ○ CUSTOM BED LENGTH ○ 60" or ○ 85" ○ JENX DREAMA POSITIONING SYSTEM ○ LOW AIR LOSS MATTRESS 	SleepSafe Beds®		
When ordering your <i>InSIGHT™ Bed</i> , please select the following choices:				
1. Finish:	O OAK O MAPLE O CHERRY O MAHOGANY O ALDER			
	SOLID COLOR _White _Blue _Green _Red _Yellow _Orange _Purple _Pink _Lt. Brown	n _Dk. Brown _Gray _Black		
2. Padding:	PADDING COLORCreamGrayLight BlueTaupeDark BluePink			

Safety First

SleepSafe® Beds are a proven, ultra-safe solution to address problems with entrapment and falls. Unique in the industry, these beds meet or exceed FDA guidelines on the **7 Zones of Entrapment**.

FOUR MODELS of SleepSafe® Beds are available in TWIN, FULL and QUEEN sizes with safety rail height protection ranging from 8¾" to 53" above the mattress depending on the configuration. Appropriate SAFETY RAIL HEIGHT is based on the activity level of the user, with more active users requiring additional height.



Resists Burrowing

The **SleepSafe**® precision-fit custom

mattress fits tightly within the bed frame.

difficult to burrow deeply enough to create

Even with aggressive compression, it is

Premium Memory Foam Mattress

SleepSafe® Beds use CertiPUR-US® compliant foam products for our mattresses featuring gel technology. They provide optimized pressure redistribution, rapid heat dissipation, and increases in breathability compared to standard foam. Mattress coverings are fluid resistant and support up to 350 pounds.

New!

Easily Moved

and can be locked.

Casters provide mobility

SuperSceen Top enclosure is now standard for Extension Beds.



Safest Yet!

The SleepSafe® Extension Bed in bed safety height another 1634" ab safety rail height of the SleepSafer' The bed may be ordered with or wi a top enclosure. A retrofit kit for mo existing SleepSafer® Beds is availal



Easy-access Safety Rails

Full-length wood safety rails are both hinged and removable for quick and easy accessibility.

Choose Clear or Mesh Windows

Clear, break-resistant windows, or our new mesh windows, provide visibility for both the occupant and caregiver. The need to see and be seen is important for comfort and safety.

TRANSFER HEIGHTS

Please consider the TRANSFER HEIGHT before ordering your bed. It is the distance from the top of the side panel (below the safety rail) to the floor. SleepSafe® - LOW - 20¾" or 22¼" Transfer Height SleepSafe® II - MEDIUM - 30½" or 32" Transfer Height SleepSafer® - TALL - 24½" or 26" Transfer Height

New Soft-Sided Solution!

InSIGHT[™] Bed by SleepSafe® Beds

Offers a 360° see-through screen surround and top for full visibility. Includes a water resistant Mattress/Bunkie Board combo, tubing channel, casters, anti-tip hardware, and padded corners.



SleepSafe® Bed Model Measurements

Use this chart to compare the different measurement specifications from model to model.

SleepSafe Bed, LLC 3629 Reed Creek Dr. Bassett, VA 24055


Toll Free: 866.852.2337 SleepSafeBed.com

SleepSafe® Bed	Low Bed	Medium Bed	Tall Bed / Extension Bed	
Specifications	SleepSafe®	SleepSafe® II	SleepSafer®	
FOUNDATION SAFETY RAIL HEIGHTS The distance from	the top of mattress to the top o	f the safety rail.		
Fixed w/Box Spring (Frame attachment choice)	12" or 8¾"	17¾", 14" or 10"	31" or 24¾" / 47¾" or 41½"	
Fixed w/Bunkie Board (Frame attachment choice)	171/4" or 14"	27" or 21¾"	36¼" or 30" / <mark>53" or 46¾"</mark>	
Manual Articulation (Frame attachment choice)	15" or 11¾"	21", 17" or 13¼"	34" or 27¾" / 50¾" or 44½"	
Manual Electric Combo (LO to HI Settings based on frame attachment)	LO 11½" to EVEN with rail	LO 21¼" to HI 9½" or LO 16" to HI 4¼"	LO 30½" to HI 18¾" / LO 47½" to HI 35¾"	
Manual Hi-Lo (LO to HI Settings based on frame attachment)	LO 11½" to EVEN with rail	LO 21¼" to HI 9½" or LO 16" to HI 4¼"	LO 30½" to HI 18¾" / LO 47½" to HI 35¾"	
Manual Hi-Lo Combo (LO to HI Settings based on frame attachment)	LO 11½" to EVEN with rail	LO 21¼" to HI 9½" or LO 16" to HI 4¼"	LO 30½" to HI 18¾" / LO 47½" to HI 35¾"	
Electric Articulation (Frame attachment choice)	15" or 11¾"	21", 17" or 13¼"	34" or 27¾" / 50¾ or 44½"	
Electric Hi-Lo Combo (LO to HI Settings based on frame attachment)	LO 11½" to HI 1¼" below rail	LO 21¼" to HI 11" or LO 16" to HI 5¾"	LO 30½" to HI 20¾" / LO 47½" to HI 37¼"	
FOUNDATION MATTRESS TO FLOOR HEIGHTS The dis	stance from the top of the mattr	ess to the floor.		
Fixed w/Box Spring (Frame attachment choice)	23¾" or 27"	28", 32" or 35½"	23¾" or 30"	
Fixed w/Bunkie Board (Frame attachment choice)	181/2" or 213/4"	18¾" or 24"	181/2" or 243/4"	
Manual Articulation (Frame attachment choice)	20¾" or 24"	25", 28¾" or 32½"	20¾" or 27"	
Manual Electric Combo (LO to HI Settings based on frame attachment)	LO 24¼" to HI 36"	LO 241/2" to HI 361/4" or LO 293/4" to HI 411/2"	LO 241/4" to HI 36"	
Manual Hi-Lo (LO to HI Settings based on frame attachment)	LO 24¼" to HI 36"	LO 241/2" to HI 361/4" or LO 293/4" to HI 411/2"	LO 241/4" to HI 36"	
Manual Hi-Lo Combo (LO to HI Settings based on frame attachment)	LO 24¼" to HI 36"	LO 241/2" to HI 361/4" or LO 293/4" to HI 411/2"	LO 241/4" to HI 36"	
Electric Articulation (Frame attachment choice)	20¾" or 24"	25", 28¾ or 32½"	20¾" or 27"	
Electric Hi-Lo Combo (LO to HI Settings based on frame attachment)	LO 24¼" to HI 34½"	LO 24½" to HI 34¾" or LO 29¾" to HI 40"	LO 241/4" to HI 341/2"	
FLOOR TO RAIL HEIGHT The distance from the floor to the top of the rail.	35¾"	45¾"	54¾" / 72"	
OVERALL OUTSIDE DIMENSIONS				
Twin	42¼"W x 80½"L x 40¾"H	42¼″W x 80½″L x 50¾″H	42¼"W x 80¼"L x 57"H (72¾"H w/Ext.)	
Full	56"W x 80"L x 40¾"H	56"W x 80"L x 50¾"H	56"W x 80"L x 57"H (72¾"H w/Ext.)	
Queen	65"W x 80"L x 40¾"H	65"W x 80"L x 50¾"H	65"W x 80"L x 57"H (72¾"H w/Ext.)	
MATTRESS SIZE Width without pads / Width with pads / Leng	ŋth			
Twin / Full / Queen	38"W / 36"W / 77"L	52"W / 50"W / 77"L	61"W / 60"W / 77"L	
UNDER BED FLOOR CLEARANCE				
Standard Models	8"	8"	7¾"	
TRANSFER HEIGHT The distance from top edge of the folded down safety rail to the floor.	22¼" with Safety Rail On 20¾" with Safety Rail Off	32" with Safety Rail On 30½" with Safety Rail Off	26" with Safety Rail On 241/2" with Safety Rail Off	







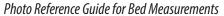



SleepSafe®

SleepSafe® II

SleepSafer®



needed. A bunkie board allows 6" of addtional safety rail height.

ARTICULATED — The mattress rests on an ELECTRIC or MANUAL frame that has both head and knee articulation, operated by a hand-held remote or crank.

HI-LO — The mattress rests on an ELECTRIC frame that offers mattress height adjustability as well as head and knee articulation. The surface of the mattress can be adjusted from 22¾" to 38¾" from the floor.

Accessories

PADDING

Custom-fit padding provides a layer of protection from the bed frame. Made with 3/4" high-density foam, it is covered with water-resistant nylon. Pads can either surround the windows or cover the windows.

IMPORTANT: When required, padding needs to be indicated when placing the order for your bed.

DREAMA® POSITIONING SYSTEM

A 24-hour postural management support system that can be used in a SleepSafe bed, and in most cases, be retrofitted to an existing bed.

LOW AIR LOSS MATTRESS

An Alternating Air-pressure Mattress from Blue Chip Medical, can now be fitted for use in a SleepSafe® Bed.

DUAL VIEW SAFETY RAILS

May be ordered for a SleepSafer Tall Bed with safety rails on **both** sides for dual access and visibility.

IV KIT / TUBING CHANNEL

An IV pole kit and a 2"x18" channel for medical tubing, may be ordered with your bed.



Finish Choices

0

SleepSafe® Beds can be finished in wood or color choices:



